




INDEPENDENT X-RAY  
QUALITY ASSURANCE

Use with

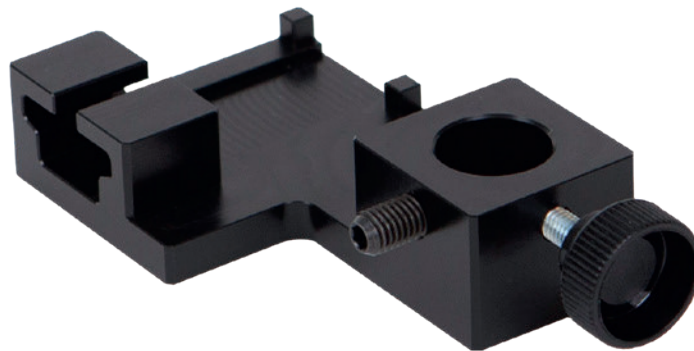

## Philips Mako Holder | Specifications

# Simple Positioning of Mako R/F Probe on Philips X-ray Systems

The RTI Philips Mako Holder allows for easy and precise positioning of the Mako R/F Probe into the correct reference position for dose measurement on Philips Interventional/Angio/Cath lab X-ray systems.

The holder fits to the Philips EDL Tool, which has a cylindrical protrusion to which the holder attaches. The Mako Probe is connected to the Mako Base Unit via USB cable.

### Application

This Philips EDL Tool for Interventional X-ray systems protects the detector from X-ray during QA procedures. The RTI Philips Mako Holder can be attached to the Philips EDL Tool on the cylindrical protrusion, with a locking mechanism to hold it in place. The Mako R/F Probe fits into the holder into the correct position for dose measurement.

The Philips EDL Tool comes with a white holder. To use Mako with the Philips EDL Tool, simply replace the white tool with the RTI Philips Mako Holder.

### Specifications

Article number:	9765033-00
Article name:	Philips Mako Holder
Dimensions:	122 x 63 x 27 mm
Weight:	92 g
Material:	POM
Compatible meters:	Mako Probes
Compatible systems:	Philips X-ray systems (via Philips EDL Tool)

